



# User Instruction Manual

## Hoyer® Standard Deluxe Sling

To avoid injury, read user manual prior to use.

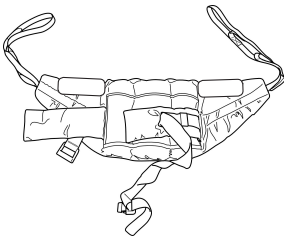

**Deluxe Standing Sling**

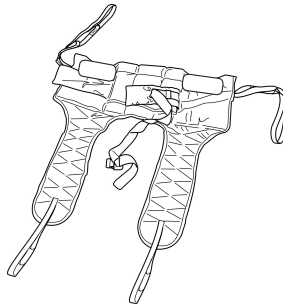

**Deluxe Transport Sling**


### **⚠ WARNING**

- **JOERNS RECOMMENDS THE USE OF GENUINE HOYER SLINGS.** Hoyer brand slings are approved by Joerns for use on Joerns brand lifts. Slings used on Joerns brand lifts, including Hoyer lifts, must be approved by Joerns.
- For the safety of the patient and carer; before using a sling a full risk assessment must be conducted to ensure that the correct sling choice, method of positioning in the sling and procedure for transfer has been determined for the patient.
- **CHECK** sling and stitching before each use. Using bleached, torn, cut, frayed or broken slings is unsafe and could result in serious injury or death to the patient.
- **DO NOT** alter slings. Destroy and discard worn slings.
- **NEVER** leave a patient unattended.
- **DO NOT** exceed the rated capacity of the sling or lift.
- **DO NOT** attempt to re-position a patient by pulling on the sling loops.


## Introduction

The Deluxe Standing sling is an easy fit sling suitable for standing and supporting. It is not a general purpose sling and clients should be assessed as to their suitability. The Deluxe Transport sling is an easy fit general purpose sling giving full support to most clients. It will not suit double amputees or lifting from a horizontal position.

## How to Use the Deluxe Standing Sling

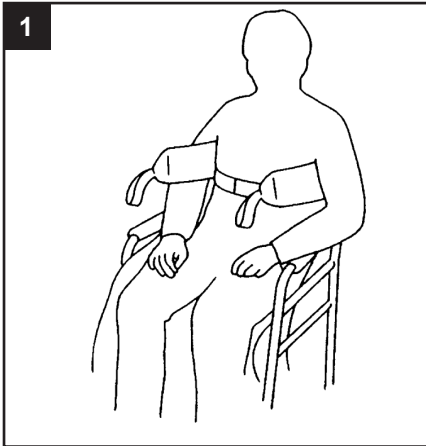

Place sling down behind client's back between base of shoulder blades and bottom of rib cage. Attach strap around the front of the client. This strap need not be tight. It is there to hold the sling in position whilst attaching the sling to the Standaid. Ensure that the client's arms are outside the sling.

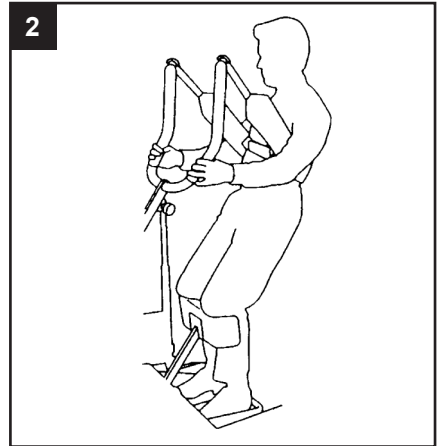

Position Standaid so that the client's feet fit in the foot trays and knees rest against the knee pad. Initially attach the longer loops of the sling to the white hooks at the end of the hand grips. The client's hands should be placed on the black hand grips or support handles. (Once the client is used to the lift, the shorter loops may be used to achieve a more upright position). Then push the up button on hand control to raise the client to the standing position.

**NOTE:** THIS SLING IS NOT A GENERAL-PURPOSE SLING AND WILL NOT BE SUITABLE FOR ALL CLIENTS.


### **TIPS FOR USE:**

CLIENT SHOULD BE ASKED TO LEAN BACK INTO THE SLING FOR THE DURATION OF THE LIFT.


TALLER CLIENTS MAY NEED TO USE SHORTER LOOPS ON THE SLING FROM THE START.

SHOULD THE CLIENT START TO LOSE TRUNK CONTROL YOU WILL FIND THAT THE SLING WILL SLIP UP BENEATH THE ARMPITS MAKING THE LIFT UNSAFE. IN THIS SITUATION THE TRANSPORT SLING CAN BE USED PROVIDING A RISK ASSESSMENT IS COMPLETED.


## How to Use the Deluxe Transport Sling

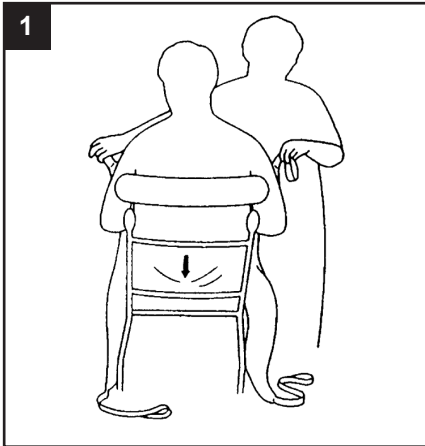


Feed the sling down the back of the client leaving the top of the commode aperture at the base of the spine. Check the sling is square across the client's base.

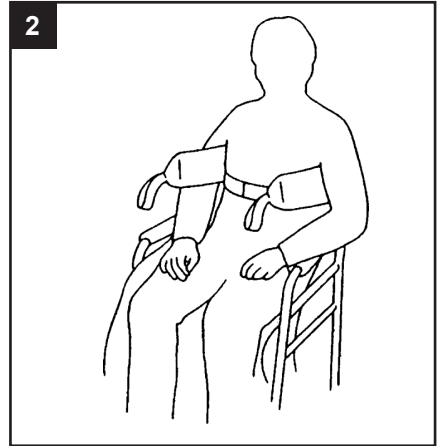

Attach the strap around the chest of the client. This strap need not be tight, it is there to hold the sling in position whilst attaching the sling to the Standaid. Ensure client's arms are on the outside of the sling.

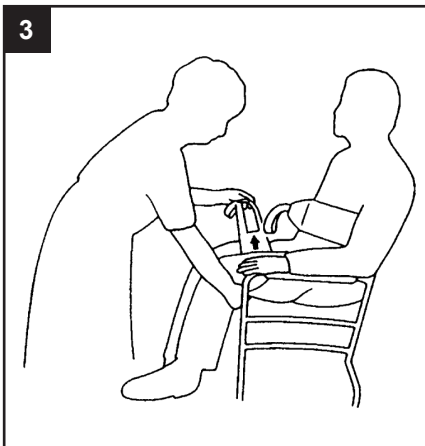

Raise the client's leg and feed the leg strap under, and then up between the legs. Ensure that the sling is not twisted or creased under the thigh. Repeat this procedure for the other leg. Attach the loops on the top of the sling to the white hooks at end of hand grips.

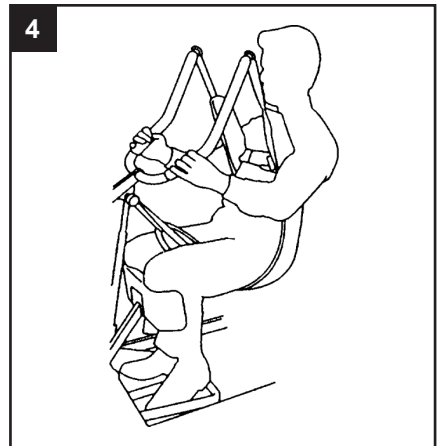

Without crossing loops on the leg straps, attach to left and right black knobs on the lifting fork(s). You would normally use the middle of the three loops. You are then ready to raise the client from the seat and transport. Once clear of the seat, the client will be more comfortable if lowered back to sitting position.

**ENSURE THAT YOUR CLIENT'S HEAD IS SUPPORTED AT ALL TIMES**

## Technical Specifications

STANDARDS	EN ISO 10535, CE
SAFE WORKING LOAD	440lbs / 200kg

## Sizing & Part Numbers

SIZE	STANDING	TRANSPORT
PEDIATRIC (BROWN)	SPECIAL	SPECIAL
SMALL (RED)	NA13499N	NA13502N
MEDIUM (YELLOW)	NA13500N	NA13501N
LARGE (GREEN)	NA13498N	NA13503N
EXTRA LARGE (BLUE)	SPECIAL	SPECIAL

## Washing Instructions



Machine wash at 185°F.



DO NOT wash with bleach. Bleach will damage the sling's material and make it unsafe for use.



Cool tumble dry, air dry or dry at very low temperature.



DO NOT dry clean.

### **⚠ WARNING**

**Slings can suffer damage during washing and drying and should be checked carefully before each use.**

Additional slings are available to suit individual needs. Contact your Joerns representative for more information.