



Optional Liners
shown

 special
Tomato
JoggerTM

ITEM# 62000000
REVISION: 01/14

UMJOG01-14LIT

In the interest of your priceless cargo, please read these warnings before using your Tomato Jogger.

 **WARNING: AVOID SERIOUS INJURY FROM FALLING OR SLIDING OUT. ALWAYS USE THE SEAT HARNESS.**

NEVER LEAVE CHILD UNATTENDED.

IMPORTANT: KEEP THESE INSTRUCTIONS FOR FUTURE REFERENCE. A CHILD'S SAFETY IS YOUR RESPONSIBILITY. THIS JOGGER IS SUITABLE FOR ONE CHILD. ALWAYS USE THE ATTENDANT STRAP. TOTAL WEIGHT LOAD RECOMMENDED FOR THIS JOGGER IS 110LBS OR 50 KG.

GENERAL INFORMATION & SAFETY

- Front wheel must be kept in-line for exercise.
- Follow all instructions carefully. Incorrect use will cause damage to push chair.
- This product has been designed to take the maximum weight of 110lbs (50kg).
- Children should be harnessed in at all times and should never be left unattended.
- **IMPORTANT - It is important for your child's safety that the harness be correctly fitted and adjusted. If not adjusted correctly in accordance with the instructions, the stability of the Jogger will be compromised.**
- The harness and seat belt are not a substitute for proper adult supervision.
- The child should be clear of moving parts while making adjustments.
IMPORTANT - DANGER: It is important that care is taken to ensure that children are kept clear of the Jogger when folding or unfolding. Pinch points and scissoring actions are unavoidable during these operations.
- Do not carry additional children, goods or accessories in or on this push chair except as permitted in this manual. They may cause the Jogger to become unstable or place unneeded stress on the Jogger that could lead to breakage.
- Overloading, incorrect folding and the use of accessories, e.g. child seats, bag hooks, rain covers, buggy boards, etc., other than those approved by the manufacturer may damage or break this Jogger.
- Any damage caused by the use of accessories not supplied by manufacturer will not be covered by the terms of our warranty.
- Never leave child in the Jogger when ascending or descending stairs or escalators or when travelling on other forms of transport.
- The Jogger should always open and fold easily. Do not force it. It may be necessary to re-read the instructions.
- Negotiate curbs and rough ground carefully. Repeated impact will cause damage to the push chair.
- **Never allow children to stand on the footrest.**
- Do not leave the Jogger exposed to heat, i.e., by a radiator or in direct sunlight.
- Do not place items on the top of the hood as it may cause damage to the canopy.

MAINTENANCE AND REPAIR

This Jogger requires regular maintenance by the user. Before each use check all rivets and connecting devices for tightness and security. Inspect all brakes, wheels, and tires and replace or repair if necessary. Check all safety devices for correct operation. These should be free to move at all times. Never continue to use a product that is not structurally sound.

- We recommend a service every 12 months.
- If wheels squeak, lubricate sparingly with Teflon or silicone coating. Do NOT use oil or grease based products, as this will attract dirt, which will clog the movement.
- Only manufacturer replacement parts should be used. It may be unsafe to use parts not supplied by manufacturer.

CLEANING

- Clean frame with a damp cloth and a mild detergent and dry.
- If the Jogger parts have been exposed to salt water, we recommend that they be rinsed down with fresh (tap) water as soon as possible.
- The canopy may be sponged lightly using a damp cloth and a mild detergent. The seat cover and head cushion may be washed. Line dry before reuse.
- Do not fold or store the product while wet and never store in a damp environment, this can cause mildew to form.

The manufacturer reserves the right to make design changes to any of its products as part of its continuous design improvement program.

INSTRUCTIONS

IMPORTANT:

- READ THESE INSTRUCTIONS CAREFULLY BEFORE USE AND KEEP THEM FOR FUTURE REFERENCE.
- YOUR CHILDREN'S SAFETY MAY BE AFFECTED IF YOU DO NOT FOLLOW THESE INSTRUCTIONS.
- APPLY BRAKE WHEN PARKING THE JOGGER.


Contents

2	Frame
3	Front Wheel Attachment
3	Front Wheel In-line Locking
4	Rear Wheel Attachment
4,5	Safety Harness
6	Parking Brake
6	Recline Seat
7	Removing Seat Fabric
7,8	Shopping Basket
8	Canopy Attachment
8	Handle Height Adjustment
9	Folding & Optional Liners

Frame

Unpack Jogger from box and set wheels aside.

1. Push foot plate of Jogger away from handle by pushing curved crossbar up and out.

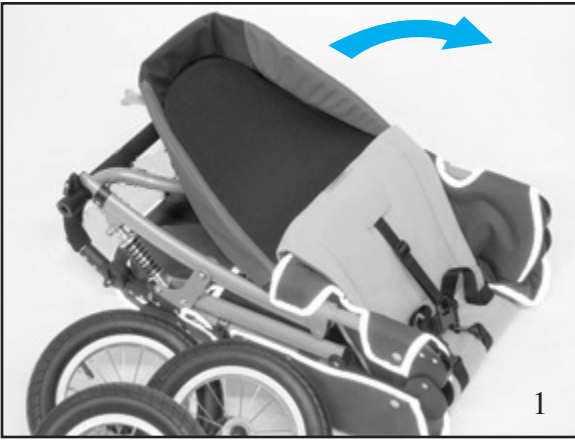


2. Audible click will be heard when Jogger frame is completely unfolded.




3. Flip Jogger to upright position.



4. Clip lower belting around frame to secure seat fabric.

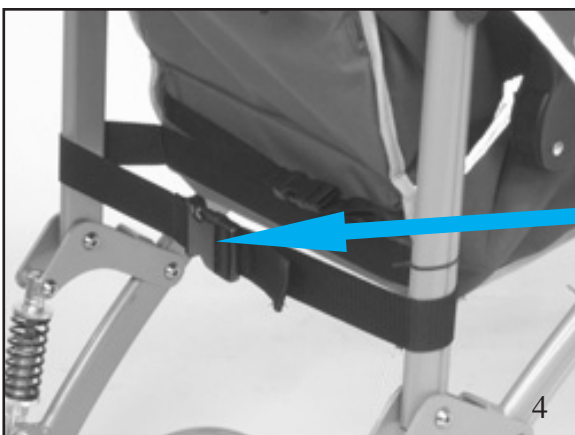

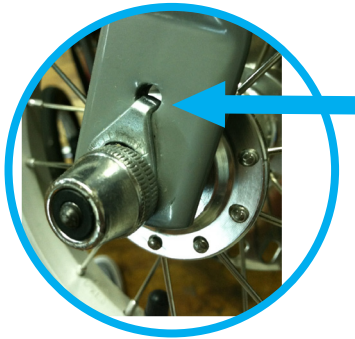


5

*Safety Lock - both sides

Front Wheel

To install front wheel, lay jogger on it's back so handle and rear axles are touching the floor.

5. Engage safety lock on both sides of the fork as shown in image 5. Check before each use.

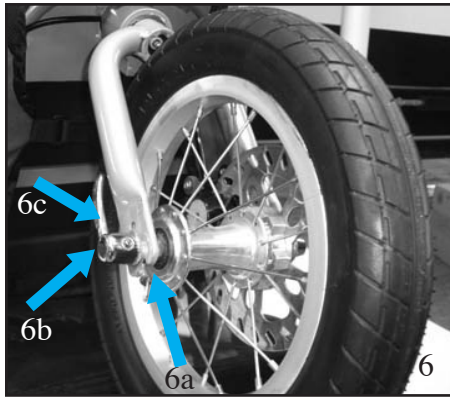

6a. Insert the Front Wheel into the fork tips by aligning the axle flats with the fork tip flat surfaces.

6b. Tighten by turning the wheel lever and securing the other side.

6c. Close the wheel lever making sure that it is flush with the tubing. Be sure not to over tighten.



7

Front Wheel In-line Locking



Front wheel must be kept in-line for exercise


7. The in-line locking mechanism is located behind the front wheel pivot point. **Pull and rotate knob to lock and unlock.** Image shows unlocked.



8

8. Image now shows the wheel locked in-line. Pull and rotate mechanism to unlock.

Rear Wheels

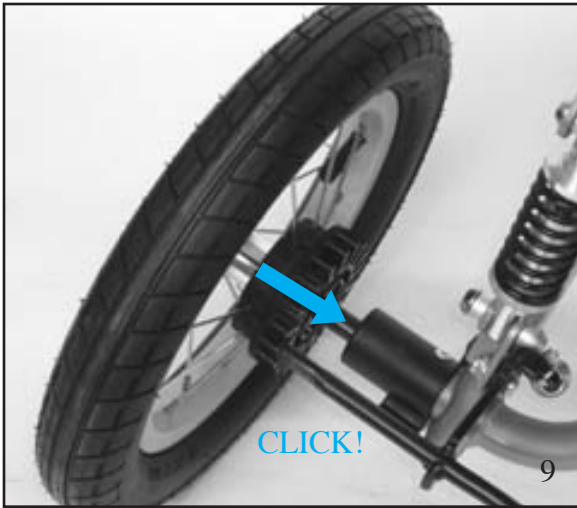

9. To install, gently slide the rear wheel's axle into the frame's supporting axle hole until you hear it click into place.




10. To remove, simply push down the spring lock tab while pulling out on the wheel.




Safety Harness

Adjusting safety harness:

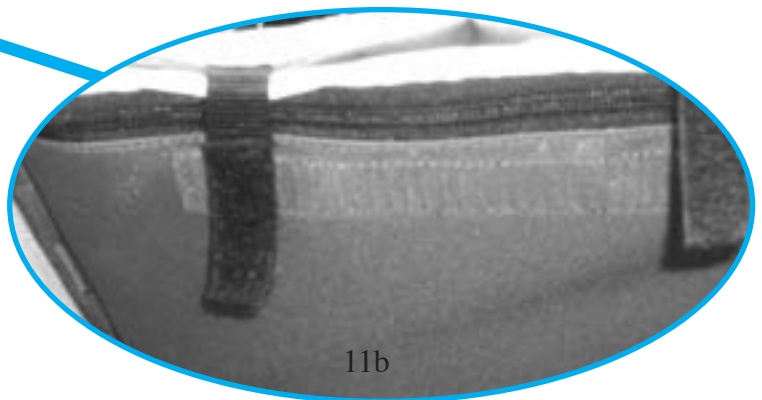


11a. Unclip hooks.

11b. Remove fleece lining (Velcro tabs shown above).

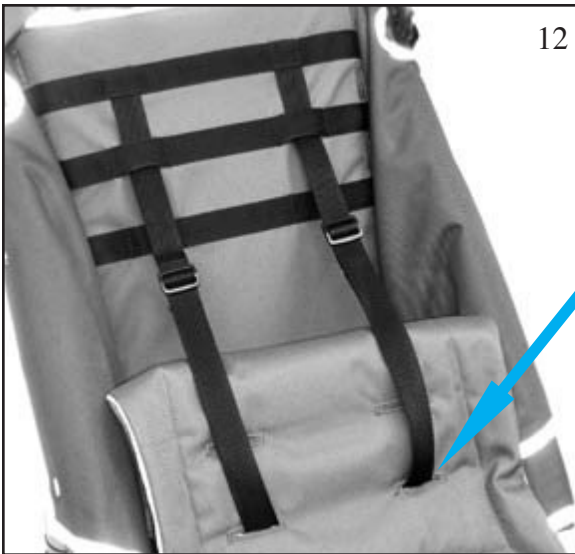


12. Unthread harness straps from fleece lining.

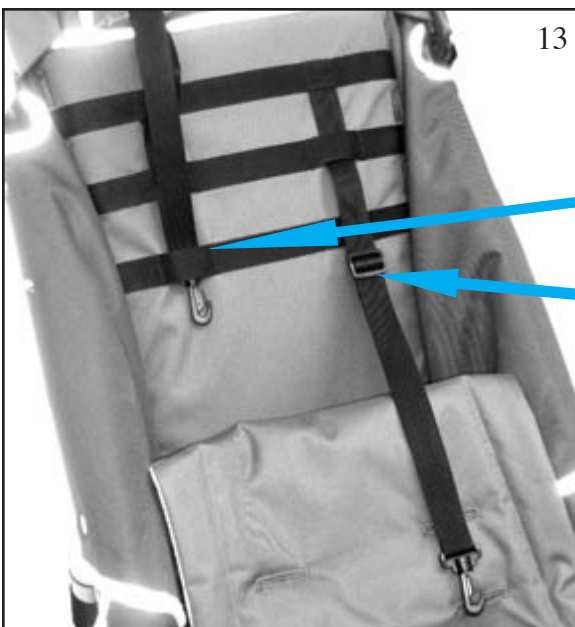

13a. Adjust shoulder straps to desired height using the loops on the seat back.

13b. Slide the plastic guide up and down to adjust strap length.

Once fitting is complete, reassemble the fleece lining.

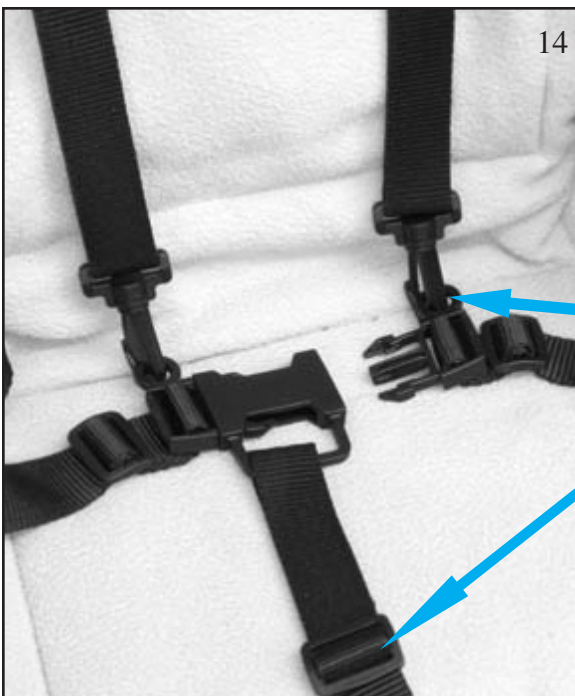


14a. Clip both shoulder straps to the eye holes located on the buckle mechanism.

14b. To adjust the length of the center strap, move the plastic guide up and down.

Rear Wheel Parking Brake


To apply the parking brake:

15. Press down on the foot bar.

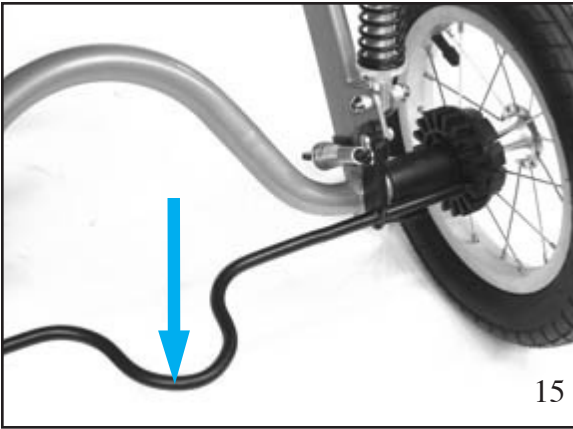

Disengage the brake:

16. Pull up on the foot bar.

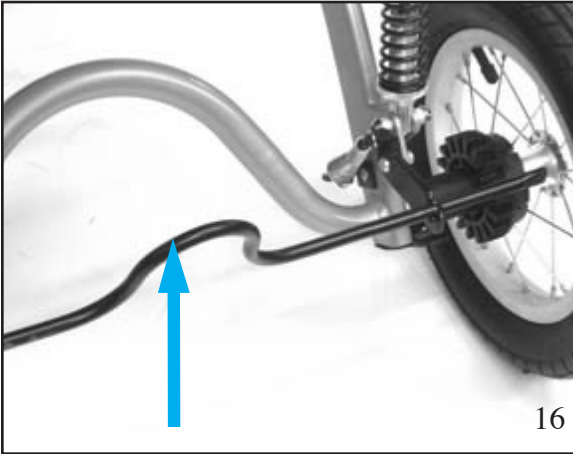

Reclining Seat


17. Simply adjust the recline strap length to desired degree of recline.




18. Unlock the buckle shown in Figure 17 for maximum recline.



Removing Seat Fabric

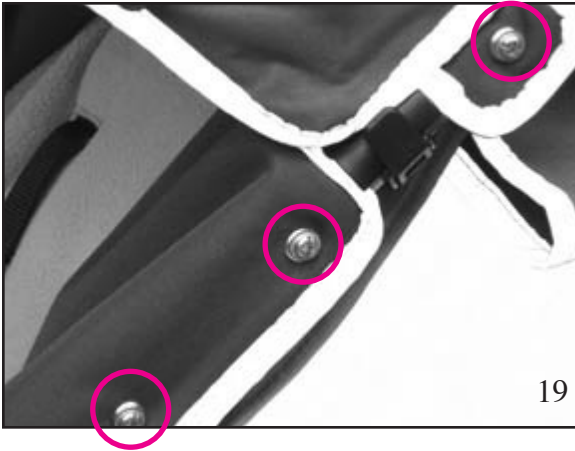


19. Using a phillips head screw driver unscrew hardware from the frame.

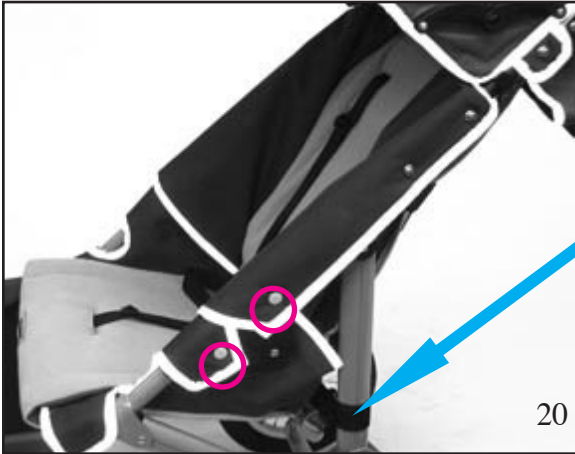

20. Unsnap fasteners on both sides of jogger.


20a. Unlock seat strap from frame, buckle is suspended across the back of the jogger.

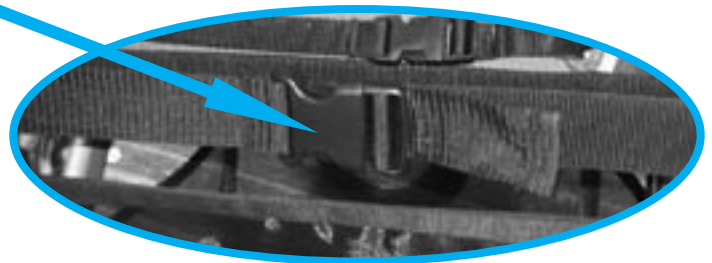

21. Unloop Velcro Strap through the D-ring, underneath the foot well.


NOTE: REMOVE ALL RIGID PARTS BEFORE WASHING
Joggers are shipped with Seat pre-assembled. Assembly/ Removal instructions are given if any pre-assembly items are removed for cleaning or replacement. "SEE PAGE 10"

Shopping Basket

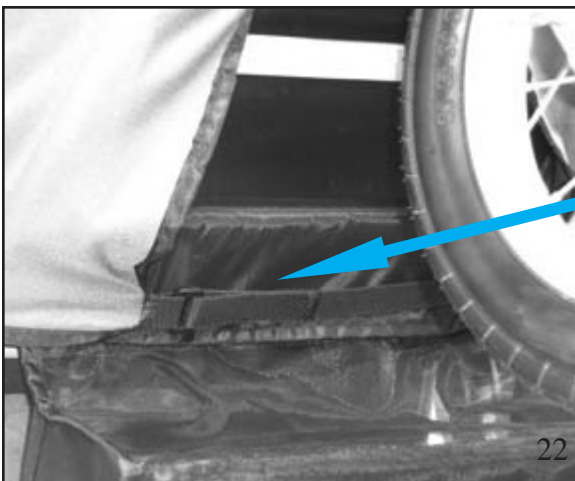

22. Secure front basket flap by snapping around the lower foot well crossbar.



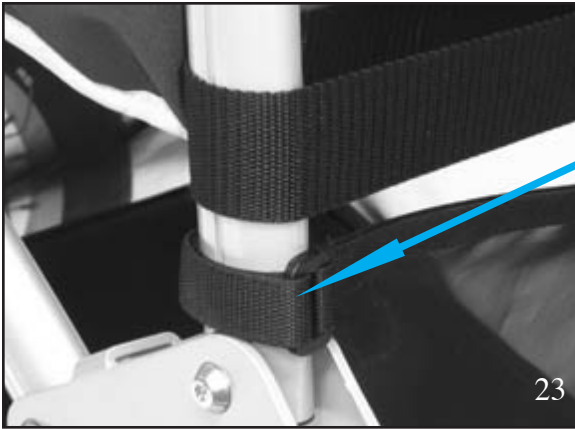


23. Wrap rear basket straps around the frame and loop back through the D-ring. Secure the Velcro strap.

Note: The D-ring loops are at the rear of the basket.

Canopy




24. Fasten the canopy to the portion of the upper frame that is not covered by fabric. Small holes in the crest of the frame tubing accept the canopy. Be sure there are NO obstructions between the canopy clip and the frame. There is a small plastic nipple on the canopy mechanism that fits into the hole in the frame.

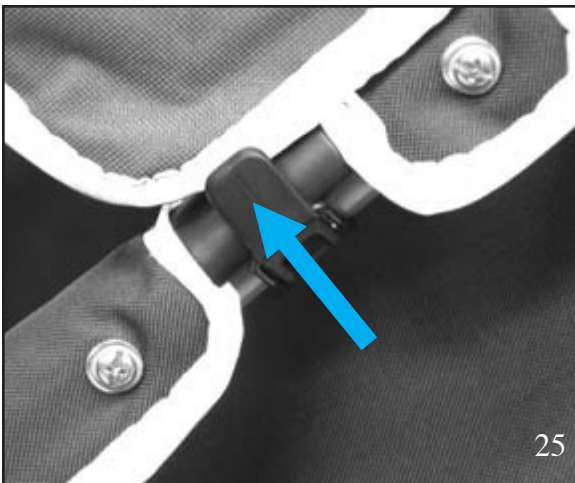

25. Hook the wire on the back side of the canopy mechanism and then close the lever by pressing it to the frame (Figure 25). **DO NOT FORCE TO CLOSE!** Double check for obstructions and be sure nipple is inserted into frame.

Handle Height

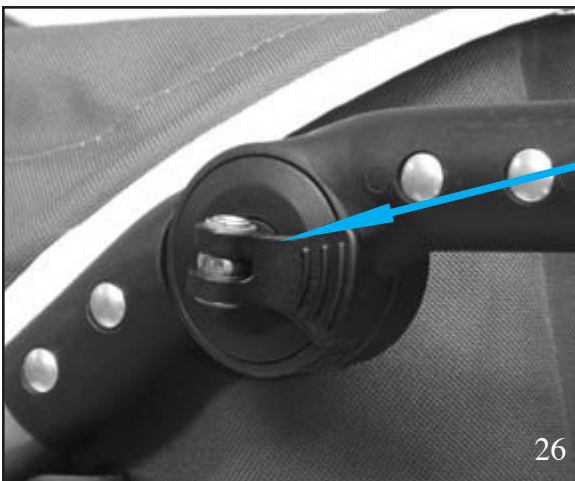

26. Unlock lever on both sides of handle. Adjust to height. Then relock.





Folding

To ease Jogger Folding:

- Fully recline seat back
- Disengage Velcro Between Canopy & Seat Back
- Push handle all the way up
- Unlock front wheel swivel lock
- Undo 5 point harness



27. Squeeze mechanism located below the handle. (Mechanism located on Both Sides)

Fold jogger back.



Recline strap secures Jogger in folded position.

Optional Special Tomato Liners



28. Back Liner:

Using the provided straps that come with the back liner – Feed straps through the slots in the jogger seat back. Extra slots provide adjustment. (Do not feed Jogger shoulder straps through liner.)

29. Seat Liner:

Use attachment straps provided with the seat liner to fasten the optional back pad to the jogger.

Then feed the original jogger seat buckles up through the seat liner. [White Buckles]

Attendant Strap

- Make sure safety wrist strap is secured properly around your wrist whenever your child is in the jogger.

Washing

- You can machine wash your Jogger seat in cold water with a mild detergent. To minimize shrinkage, simply reattach the seat to the frame to dry. (DO NOT put seat in the dryer)
- Be sure to remove the rigid back board from the seat before washing. Unzip seat back and remove rigid back support.
- You may clean the Jogger frame using a soft, dry cloth.



LIMITED WARRANTY

Your push chair is guaranteed to be free from any manufacturing defects for a period of 2 years from the date of purchase under normal use and provided it is used in compliance with the operating instructions. This warranty extends only to the original retail purchaser and is only valid when supplied with proof of purchase. Please retain proof of purchase with this Limited Warranty.

Should a repair be needed please contact the retailer from whom the product was purchased.

Please note that the warranty will not be applicable if:

A fault is caused by misuse or poor maintenance.

Repairs are carried out by a third party.

The Jogger is faulty due to general wear and tear which is the result of everyday use.

The Jogger is damaged as a result of an accident.

The manufacturer is not satisfied that the stated warranty terms and conditions have been met.

AUTHORISED EUROPEAN REPRESENTATIVE

GLOBAL CAREHAB
ROSENVANGS ALLE 32
8260 VIBY J
DENMARK
+45 24 65 55 50



WARNING: Prevent serious injury or death from falls or sliding out

- Always secure child in the restraint
- Never leave child unattended
- Always keep child in view
- Max User Weight 110lbs or 50kg



Bergeron By Design
15 South Second St.
Dolgeville, NY 13329
1.866.529.8407



T | 315 429 8407
F | 315 429 8862

© 2014 Bergeron By Design

The information contained herein is subject to change without notice. Bergeron By Design shall not be liable for technical errors or omissions herein.