

PR514-MLA

MaxiComfort with Twilight Positioning PR514-MLA

> Discover weightlessness in the MaxiComfort Cloud with Twilight positioning. Our patented threemotor technology offers an array of rejuvenating options to cradle the body and ease the mind. Achieve extreme Zero Gravity and TV positions while enjoying the world's finest rejuvenating seating experience!
















* upcharge applies

PR514-MLA

| Model Number | PR-514-MLA |
|---|--------------|
| Size | Medium/Large |
| Weight Capacity | 375 LB |
| Recline Positions | MaxiComfort |
| Standard Fabrics | 8 |
| Fabric Options | All |
| Chaise Seat | Yes |
| Coil Springs | No |
| Back Type | Biscuit |
| Head Rest Cover | Yes |
| Arm Covers | No |
| Battery Back Up | Yes |
| Overall Width | 39.5" |
| | 46.5" |
| Overall Height | |
| Width Between Arms | 20" |
| Floor to Top of Seat | 20.5" |
| Seat Depth | 20" |
| Seat to Top of Back | 32.5" |
| Seat to Top of Arm | 7.5" |
| Seat Angle Raised | 15° |
| Floor to Top of Foot Rest Reclined | 25" |
| Floor to Top of Back Reclined | 27" |
| Back Angle Closed | 105° |
| Back Angle Reclined | 165° |
| Foot Rest Angle Closed | 90° |
| Foot Rest Angle Opened | 180° |
| Distance Required from Wall Reclined | 24" |
| Extension Reclined from Back of Seat to End of Foot Rest | 43" |
| Life Time Warranty, Chair Frame, Lift Frame, Recline Mechanism | Yes |
| Warranty on Electrical parts | 3 YR |
| Pro Rated Years 4 Thru 7 | Yes |
| Warranty on Mechanical Labor | 3 YR |
| Warranty on Electrical Labor | 1 YR |
| Removable Back Option Available | No |
| Smart Tek System Standard | Yes |
| Customizable * | No |
| Left Hand Control Option Available | Yes |
| Extra Pocket Option Available | Standard |
| Heat & Massage Option Available | No |
| Gel Memory Foam Option Available | No |
| Power Pillow Option Available | No |
| Auto Flex Option Available | No |
| Foot Rest Extension Option Available | No |
| Straight Lift Option Available | No |
| User Height (suggested) | 5'7"-6'2" |
| | |

* See Customize Chart for Dimension limits

Sizes may vary depending on fabric, filling material, upholstery and carpet thickness. All measurements taken on concrete floor using levels and metal rulers.

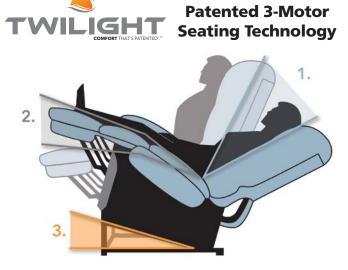

Motor #1

Adjusts back positioning

Motor #2

Controls the footrest movement and the patented movement of the rising seat box

Motor #3

Activates Twilight technology, a tilting movement that provides a weightless experience



Battery Backup System

Batteries not included. Battery Backup System is designed to bring the chair from the reclined to the seated position in the event of a power outage.

Smart Tek

Our luxurious lift and recline chairs are equipped with electronics to help assure a secure lift every time you need it.


















Literature is current at time of printing. Golden Technologies reserves the right to make changes to the product or literature at any time.

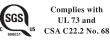



Golden Technologies • 401 Bridge Street, Old Forge PA 18518 • Tel: (800) 624-6374 • Fax: (800) 628-5165 • www.goldentech.com