



Shown with Power Elevating Seat


Standard Colors:



Crimson Red



Azure Blue

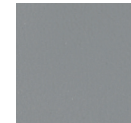

Galactic Grey

ALL-NEW COMPANION HD MODEL FOR 2024!

- 500 lb. Weight Capacity
- Extra Wide High Back Seat
- Galactic Grey, Azure Blue & Crimson Red Color Shrouds
- Tiller Mounted USB Charging Port
- 8 Amp Super Charger for Maximum Battery Life

Companion HD 3-Wheel (Heavy-Duty) GC540

GOLDEN
Companion
HD

The new Golden Companion Heavy-Duty 3-wheel scooter combines a modern design with high performance features to offer the ultimate indoor or outdoor driving experience at a weight capacity of 500 lbs. The Companion HD's rugged construction is built for stability and maneuverability with front anti-tip wheels, a 51" turning radius. It offers rear suspension for a smooth ride, and a powerful 120 Amp controller for outstanding climbing ability. It features front and rear LED lighting, an adjustable tiller with auto-lock, and a luxury high-back Captain's seat for maximum comfort when on the go, plus durable vibrant red, blue or grey color shrouds for a sporty look. With solid, non-marking tires, an operating range of 16 miles on 2-NF22 batteries, and a maximum speed of 6 mph, the Companion HD offers everything you need to enjoy the great outdoors!



Standard Features

- 500 lb. Weight Capacity
- Rear Suspension
- Powerful 120 Amp Controller
- Fully Integrated LCD Console Display that Monitors Fault Codes & Battery Energy Level
- LED Lighting Package with Front Headlights, Rear Taillights & Flashers
- Infinite Position Tiller
- USB Charging Port on Tiller
- Extra Wide High Back Seat that Swivels for Easy Access
- Large Front Basket
- 10.75" Solid Tires with Alloy Wheels
- 2-NF22 Batteries

GOLDEN Companion[™] HD

FULL-SIZE HEAVY-DUTY 3-WHEEL SCOOTER

GC540

Companion HD 3-Wheel Specifications

Model	GC540
Medicare Code	K0808
Weight Capacity	500 lbs
Drive Wheels	Rear
Maximum Speed	6 mph
Operating Range ¹	16 miles
Ground Clearance to Center Deck	6.25"
Ground Clearance to Transaxle/Motor	3.75"
Turning Radius	51"
Type Batteries / Air Travel Safe	2-NF22 / Yes
Colors	Grey, Red, Blue
Freewheel Mode	Yes
Electronic Speed Control	Yes
Electro-Mechanical Brakes	Yes
Energy Efficient Auto Turn Off (Sleep Mode)	Yes
Charger	Portable 8 Amp Off Board
Controller	120 Amp
Dimensions:	
Length	48"
Width	26"
Ground to Top of Seat	24.5" - 27.5"
Ground to Top of Deck	8"
Deck to Top of Seat	16.5" - 19.5"
Number of Seat Height Positions	4
Size of Increments	1"
Front Axle to Rear Axle	36"
Weight of Unit: (assembled)	225.5 lbs
Tires: Flat Free / Non Marking Solid	
Front	10.25"
Rear	10.75"
Front & Rear Anti-Tippers	Yes
Standard Seat: Extra Wide High Back	
Back Height (no headrest)	18"
Width x Depth	22" x 20"
Color	Black
Transportation Seat Lock	Yes, Self Locking
Seat Rotation	90 Degree either direction 360 Degrees
Optional Seat: Captain's Chair	
Back Height - no headrest	18"
Width x Depth	20" x 20"
Color	Black
Power Elevating Seat - 6" with Lockout	Yes
Warranty:	
Frame	5 Years
Drive Train	2 Years
Electronics ²	2 Years
Batteries	1 Year

- ¹ Battery range at will vary due to rider weight, drive surface, terrain, and battery type.
² Electronics warranty excludes batteries.



Digital Dashboard & User-friendly auto-tiller adjustment



Infinite adjustable tiller with ergonomic design



Front and rear LED lights and side reflectors for safety



Conveniently located drink holders



Rear suspension for a smooth ride



Extra Wide High Back Seat with Headrest for added support & comfort



Phone: 800-624-6374

Fax: 800-628-5165

www.goldentech.com

Find us on
Facebook

www.facebook.com/goldentechnologies