



If you want a great way to exercise, cross-train or just have fun, one of the **Invacare® Top End® Excelerator™ Series Handcycles** is just what you need. This upright and extremely maneuverable three-wheel handcycle is capable of speeds up to 15 mph. The smooth seven-speed Shimano® hub with its reversing drum brake delivers hands-on braking and allows the rider to tackle hills or level terrain with equal ease. The oval tube carbon steel frame construction with adjustable sliding seat makes transfers in and out easier too. There are two great models to choose from: the stocked model can be shipped right away and the custom model allows you to choose your favorite color and seat width.

The Excelerator Handcycle offers some great options too. You can order the popular mountain drive option which adds 7 extra low gears for climbing hills or the V crankset option with its wider cranks to allow more clearance for those who are tall or require additional space. No excuses...Get Cranking! – Check out our childrens' models too!

### Invacare® Top End® Excelerator<sup>™</sup> Features

Carbon steel frame and fork Adjustable seat: forwards and backwards Adjustable back height Fold over back upholstery Adjustable crank height Chain tension idler Full chain guard, reflectors and safety flag Adjustable height/angle footrest system Shimano® 7 speed hub with reversing drum brake Parking brake 175 mm crankset with one chain ring Ovalized ergonomic hand pedals Performance wheels with cruiser tires Quick release axles 3° camber Cushion and hook and loop seat restraint Weight capacity 350 lb.

# Model #: XCL/XCLSTK



- 4. Mountain drive transmission option
- 5. Frictionless ovalized handpedals
- 5. Thetiomess ovalized handpedals

#### Options

- Crank width adapter set Mountain drive transmission adds 7 low gears Quad twist shifter helps those with hand impairments Tri-pin quad handpedals QuadGrip handpedals with 3" width extensions Water bottle and cage Tow bar for wheelchair Cordless computer
- Front safety light Rear flashing light Helmet Seat restraint: Auto style fastener Mirror Crutch holder and strap Bike rack for trailer hitch Top End t-shirt

The Invacare Top End Li'l Excelerator and Invacare Top End Li'l Excelerator-2 are sized just right for children and are available in 7 speeds or 3 speeds, respectively.

### Invacare® Top End® Li'l Excelerator<sup>™</sup>-2 Features

Aluminum frame and transport ready fork Adjustable seat: forwards and backwards Adjustable back height Adjustable tension back upholstery Adjustable crank height Chain tension idler Chain guard, reflectors and safety flag Adjustable height/angle footrest system Shimano® 3 speed hub with reversing drum brake Parking brake 153 mm crankset with one chain ring Ergonomic handpedals 20" Cruiser wheels and tires Quick release axles 9° camber Cushion and hook and loop seat restraint Electric Red paint Weight capacity 175 lb.



# Invacare<sup>®</sup> Top End<sup>®</sup> Li'l Excelerator<sup>™</sup> Features

Carbon steel frame and fork Adjustable seat: forwards and backwards Adjustable back height Fold over back upholstery Adjustable crank height Chain tension idler Full chain guard, reflectors and safety flag Adjustable height/angle footrest system Shimano® 7 speed hub with reversing drum brake Parking brake 155 mm crankset with one chain ring Ergonomic handpedals

### **Options for LXC2H and LXC**

Foam covered horizontal handpedals Quad twist shifter option Mirror Water bottle and cage Helmet Front safety light Rear flashing light

20" Cruiser wheels and tires Quick release axles 9° camber Cushion and hook and loop seat restraint Choice of 11 colors Weight capacity 250 lb.

Model #: LXC

and rs 50 lb.


Crutch holder and strap Bike rack for trailer hitch Seat restraint: Auto style fastener Horizontal handpedals



## Features

1. V Crankset

The Invacare® Top End® V crankset option accommodates consumers who might have "clearance" issues. Available on the Invacare Top End Custom Excelerator Handcycle.

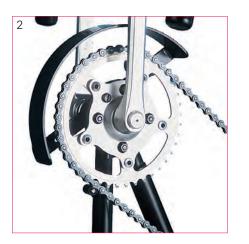

2. Mountain Drive Transmission

Available on the Excelerator handcycle Adds 7 low gears (7 low + 7 standard = 14 total) Easy shifting, especially for quadriplegics – just push in the center button Lightweight planetary transmission

A must-have for hilly terrain



### 3. QuadGrip Handpedals

QuadGrip handpedals with 3" width extensions and QuadWraps are a great option for those with limited hand function. Add the Quad twist shifter option and the mountain drive for optimal function and fun.

# Invacare Corporation

www.invacare.com

#### USA

One Invacare Way Elyria, Ohio 44035-4190 (800) 333-6900

#### Canada

570 Matheson Blvd. E., Unit 8 Mississauga, Ontario L4Z 4G4 Canada (800) 668-5324

#### Top End

www.topendwheelchair.com Sports and Recreation Products 4501 63rd Circle North Pinellas Park, Florida 33781-5914 (727) 522-8677 (800) 532-8677 Fax. (727) 522-1007

©2016

Invacare Corporation. All rights reserved. Trademarks are identified by the symbols <sup>™</sup> and <sup>®</sup>. All trademarks are owned by or licensed to Invacare Corporation unless otherwise noted. Shimano is a registered trademark of Shimano Inc. Form No. 12-509 Rev. 03/16 151174

