

## Delta Ultra-Light 1000 Full-Electric Low Bed

HCPCS: E0297



FIGURE A



FIGURE B

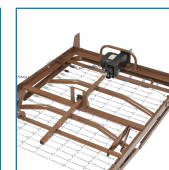


FIGURE C



FIGURE D



FIGURE E



FIGURE F

- The Delta 1000 is truly universal because the bed frame can be used with Drive's and most other manufacturers bed ends, old and new
- The transition box mounted on the foot section is key to interchangeability (Figure A)
- The transition box allows you to change the rotation of the high/low shaft
- Headboard is taller than footboard to better fit home decor (Figure B)
- The wood grain panels are break and scratch resistant
- Once installed, the high/low shaft will stay with the bed to eliminate loss or forgetting on future installations (Figure C)
- Labels and color-coded springs assure accurate installation of side rails for patient safety (Figure D)
- Unique motor is completely self-contained to reduce weight and noise
- 9 ½" to 23 ½" deck height
- Easy to set up
- Top mounted Hi/Lo crank accessibility
- Height-adjustable motor comes pre-installed on the foot section
- Caster easily moves from standard to low bed position
- The motor assembly can be installed and removed with patient in bed
- Bed, side rails and mattress safety-tested by SGS, a nationally recognized testing facility Satisfies all applicable FDA entrapment guidelines
- The motor is UL approved
- In the event of a power failure, a 9 volt battery, located in the hand control, can lower the head and foot sections nine times
- New and improved hand pendant has large, easy-to-use controls (Figure E)
- Bed ships in two cartons
- Hand crank located on the top of the footboard prevents bending for the caregiver when in the low position (Figure F)
- Channel frame construction provides superior strength and reduced weight
- Zinc-coated spring deck
- Extension kit available (model #15030EXTKIT) extends length from 80" to 84"
- Patent Pending
- Weight Capacity 450 lbs

### PRODUCTS

### Description

PRODUCTS	Description
15235	LowBed,UltraLightFull.Elec.1/c
15235BV-FR	LowBed,UltraLightFull.Elec.1/c
15235BV-HR	LowBed,UltraLightFull.Elec.1/c
15235BV-PKG	LowBed,UltraLightFull.Elec.Pkg
15235BV-PKG-1	LowBed,UltraLightFull.Elec.Pkg
15235BV-PKG-1-T	LowBed,UltraLightFull.Elec.Pkg
15235BV-PKG-2	LowBed,UltraLightFull.Elec.Pkg
15235BV-PKG-T	LowBed,UltraLightFull.Elec.Pkg

<b>SPECIFICATIONS</b>	<b>Description</b>
Bed Height Maximum With Casters	23.5"
Bed Height Minimum Without Casters	9.5"
Carton Shipping Weight	191 lbs.
Unit Dimensions	36" (W) x 88" (L)
Warranty	(Original Purchaser) 5 Years on Motor and Hand Control, Limited Lifetime on Welds and Frame, 1 Year on All Other Parts and Components
Weight	175 lbs.
Weight Capacity	450 lbs
Bed Height Min wo caster	9.5"
Bed Height Max w casters	23.5"
<b>ACCESSORIES</b>	<b>Description</b>
15030HBL	Hospital Bed Linen Kit,1 EA,10/CS, CS
15030HBC	Bedding in a Box,1 EA,5/CS, CS
<b>PARTS</b>	<b>Description</b>
15044CSET	Set of 4 Casters,Homecare Bed, 1 set ea, EA
15235	LowBed,UltraLightFull.Elec.1/c, CS
15547F-FP-1	Foot Board Panel, 15235, Delta Bed, 1/ea, EA
15547H-HP-1	HeadBoardPanel,15235/15547H,DeltaBed,1ea, EA
15044CL	Caster,w/Lockfor15544,15547,1/ea, EA
15044CNL	Caster,NoLockfor15544,15547,1/ea, EA
15541	Head Spring,15030/15033,Delta Bed,1ea, EA
15030-FRLINK	FrameLinkw/DeckScrews,15030,15005,1/ea, EA
15544HGB	Head Board Gear Box,15544,1ea, EA
15544FGB	Foot Board Gear Box,15544,1ea, EA
1503C	Locking Collar Universal Bed Shaft,1ea, EA
15544GBMB	MountingBracketForTransBox15030,1ea, EA
15033MHLNC-N	HiLoMotorConnector15033MHLN-MGBA,1ea, EA
15030UBS	Delta Bed Shaft,1ea, EA
15030MC	Cap Only for 15030 Motor,1/ea, EA
15030-MS	Motor Strap,Ultra Lite Bed,2 per PKG, PKG
15030PN-1035	HeadActuatorBarDeltaFlatTubing,1ea, EA
15030PN-1036	FootActuatorBarDeltaFlatTubing,1ea, EA
15030-7RMM	Motor Mount,7 R Serial,1/ea, EA
15030P-1006	Covers,UltraLt.Bed Okin,1/bx, BX